


# **FORNISEE**<sup>®</sup>

S Y S T E M

The only illuminated uterine manipulator for use in minimally invasive hysterectomies including LAVH, TLH, LSH, and robot-assisted GYN procedures.

**Superior visualization**  
**Precise, light-guided dissection**  
**Confident manipulation**

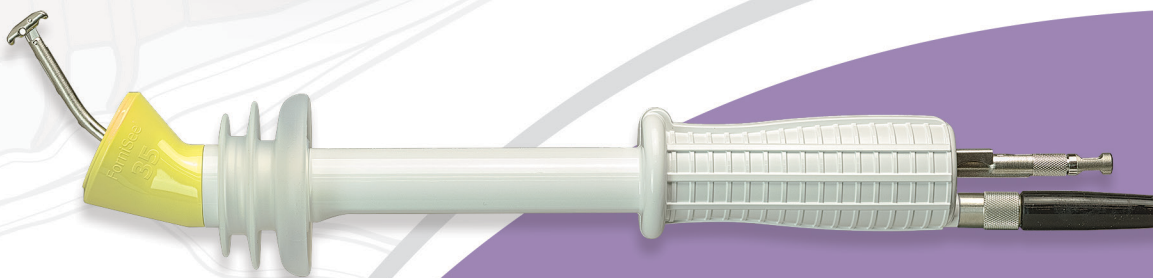

**LS SOLUTIONS**<sup>®</sup>





# FORNISEE<sup>®</sup>

## S Y S T E M

The only illuminated uterine manipulator that allows intuitive manipulation of the uterus while providing superior visualization of the fornices and precise colpotomy incisions.

Ergonomic handle allows intuitive manipulation for any size hand

### OCCLUSION

Baffled vaginal occluder maintains pneumoperitoneum

FS<sup>®</sup> LIGHT GUIDE enables precise, light-guided dissection and superior visibility

Reusable FS<sup>®</sup> SOUND with rotational anchor readily facilitates manipulation and specimen removal

Sounds: 6cm, 8cm, 10cm, 12cm

### ILLUMINATION

Illumination seal provides clear visibility of anatomical landmarks


Cervical cup fits neatly over cervix, delineates the vaginal fornices, and provides incision ridge

Cups: 30mm (green), 35mm (yellow), 40mm (blue)

### MANIPULATION


Enables anatomic tissue tensioning for uterine mobilization; Exposes broad ligaments, bladder flap, uterine blood vessels and colpotomy; Directs uterus away from ureters

#### Retroflexion

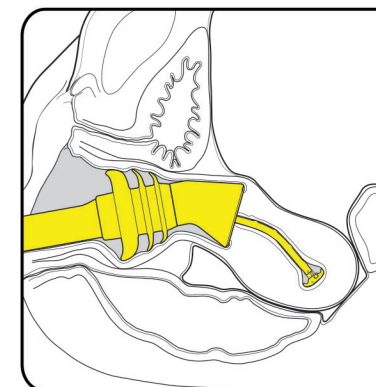


Retroflexion for anterior colpotomy and bladder flap development

#### Anteflexion

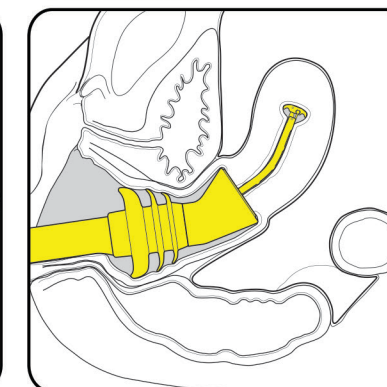

Anteflexion for posterior colpotomy dissection

#### Lateral displacement

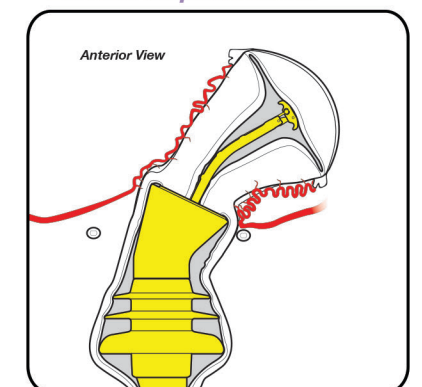

Left lateral displacement for right lateral colpotomy and right uterine artery



FS® SOUNDS			SUPPLIED: NONSTERILE	
SIZE	QTY.	REORDER	PRODUCT	DESCRIPTION
6cm	1	110056	6cm FS® SOUND - SINGLE	1 Reusable Device
8cm	1	110058	8cm FS® SOUND - SINGLE	1 Reusable Device
10cm	1	110105	10cm FS® SOUND - SINGLE	1 Reusable Device
12cm	1	110106	12cm FS® SOUND - SINGLE	1 Reusable Device
6cm, 8cm, 10cm, 12cm	1,1,1,1	110112	6cm,8cm,10cm,12cm FS® SOUNDS - SET OF 4	4 Reusable Devices

FS® DEVICES			SUPPLIED: STERILE, SINGLE PATIENT USE	
SIZE	QTY.	REORDER	PRODUCT	DESCRIPTION
30mm	6	110030	FS-30 DEVICE - (30mm Cervical Cup)	Box of 6 Single Patient Use Devices
35mm	6	110035	FS-35 DEVICE - (35mm Cervical Cup)	Box of 6 Single Patient Use Devices
40mm	6	110040	FS-40 DEVICE - (40mm Cervical Cup)	Box of 6 Single Patient Use Devices
30mm, 35mm, 40mm	2,2,2	110029	FS-30 DEVICE - (30mm Cervical Cup) FS-35 DEVICE - (35mm Cervical Cup) FS-40 DEVICE - (40mm Cervical Cup)	Combo Box of 6 Single Patient Use Devices (2 of each size)

FS® LIGHT GUIDE			SUPPLIED: NONSTERILE	
SIZE	QTY.	REORDER	PRODUCT	DESCRIPTION
7.5ft x 5mm	1	100025	FS® LIGHT GUIDE - SINGLE	1 Reusable Fiberoptic Cable

For more information, visit [www.lsisolutions.com/fornisee](http://www.lsisolutions.com/fornisee)

PN110333B 01.17.23

