


JNW URTRAC™ **INSTRUMENT GUIDE**

Dual Lumen Cannula and Obturator Kit

A multipurpose 26Fr cannula with dual
3mm and 5mm working channels.

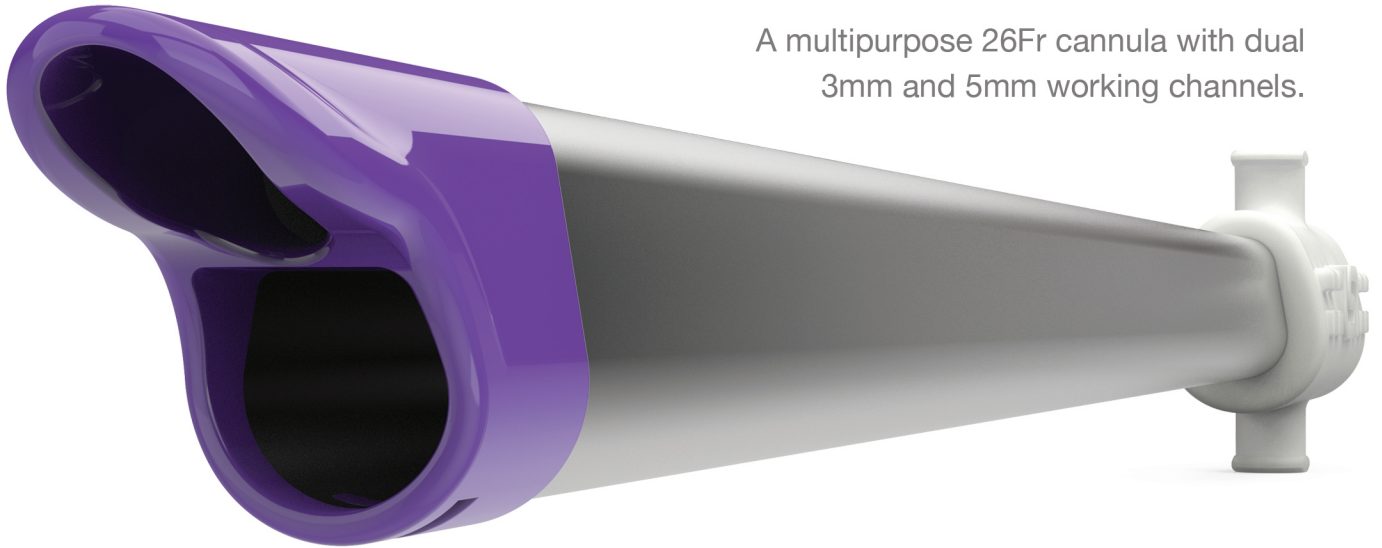

Customized Cannula Sheath
Dual Working Channels
Adjustable Flow Control


LSOLUTIONS®

JNW URTRAC™ INSTRUMENT GUIDE


Dual Lumen Cannula and Obturator Kit

A multipurpose 26Fr cannula with dual
3mm and 5mm working channels.


The JNW URTRAC™ Instrument Guide provides an atraumatic dual-lumen working channel for TUI-TMR procedures. The cannula sheath is customized to facilitate gentle insertion while providing access for simultaneous use of a 5mm shafted instrument and a 3mm flexible scope. The included dual obturator enables smooth passage of the JNW to the targeted site. The rotational seal allows for continuous flow during instrument exchanges. The JNW is compatible with the LSI SOLUTIONS® RD®, RD180®, and Ti-KNOT® devices.

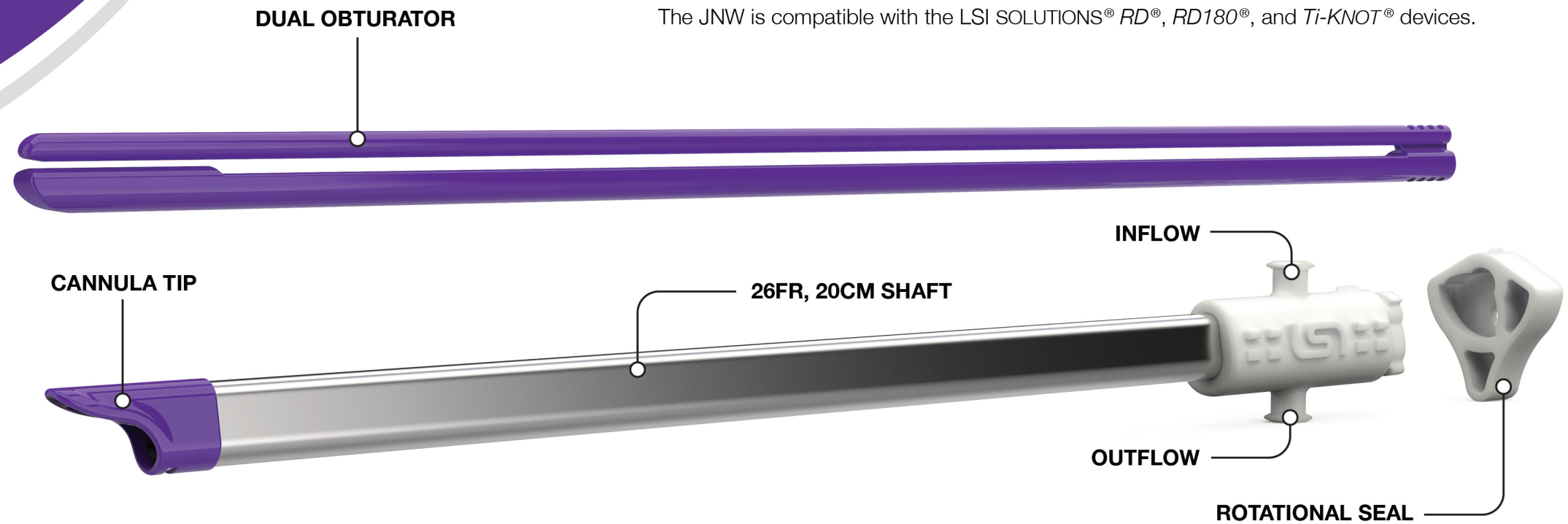


TUI-TMR: Transurethral Incision with Transverse Mucosal Realignment

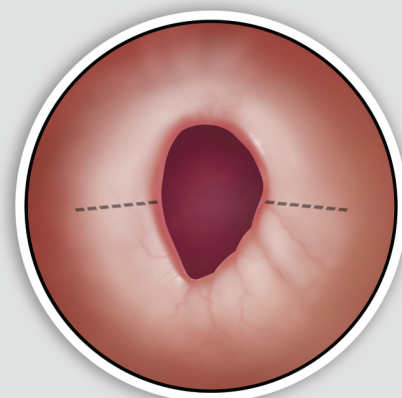


Bladder neck stricture with incisions planned at the 3 o'clock and 9 o'clock position

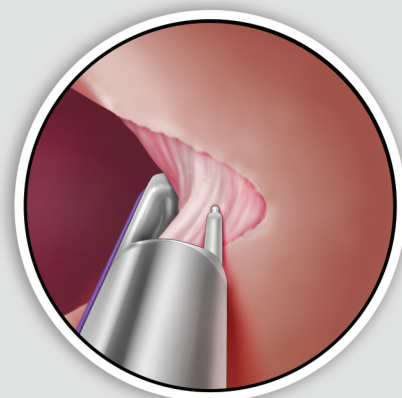


Bladder-side tissue bite taken with the RD® DEVICE at the 3 o'clock position

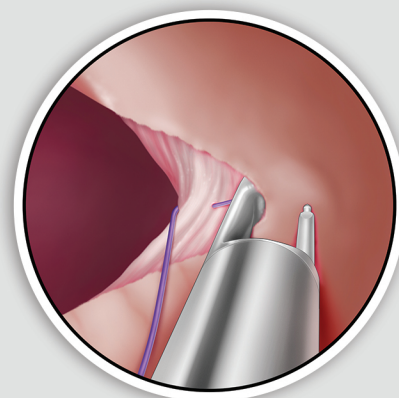


Urethral-side tissue bite taken with the RD® DEVICE at the 3 o'clock position

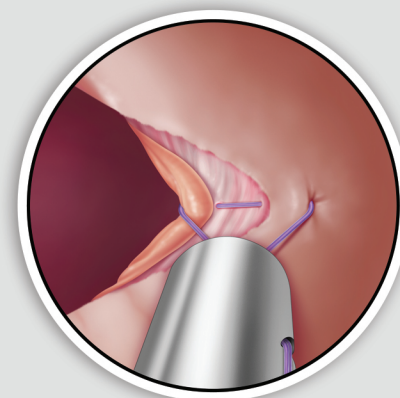

Mucosa begins to be realigned as suture is fastened with the Ti-KNOT® DEVICE

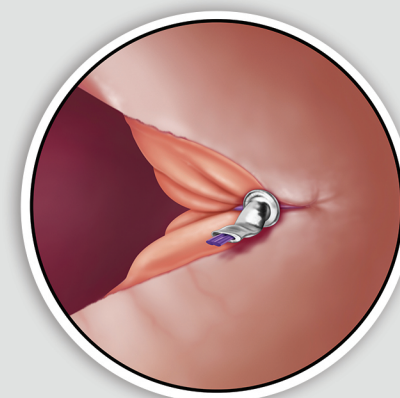


Mucosal realignment completely secured with Ti-KNOT® FASTENER


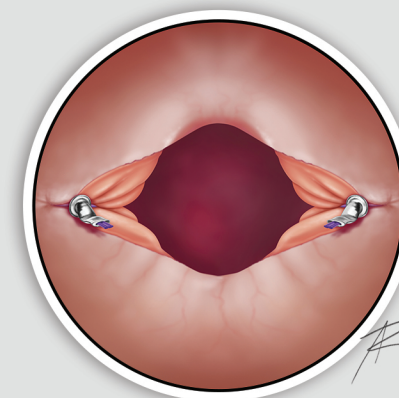



Mucosal realignment completed at both the 3 o'clock and 9 o'clock position



ACCESS

| JNW URTRAC™ PRODUCT ORDERING | | SUPPLIED: STERILE | |
|---|------------|-------------------------------------|---|
| | REORDER | PRODUCT | DESCRIPTION |
|  | REF 081300 | JNW URTRAC™ INSTRUMENT GUIDE | Box of 6 Kits (1) Cannula Sheath Per Kit (1) Obturator Per Kit (1) Rotational Seal Per Kit |

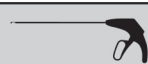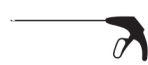SUTURING

| RD® / RD180® PRODUCT ORDERING | | SUPPLIED: STERILE | | INDIVIDUALLY PACKAGED |
|--|------------|--|---|-----------------------|
| | REORDER | PRODUCT | DESCRIPTION | |
|  x 6 | REF 021000 | 5mm RD® DEVICE | Box of 6 Devices | |
|  x 6 | REF 021100 | 5mm RD180® DEVICE | Box of 6 Devices | |
|  x 12 | REF 021030 | MONOGLIDE® RD® QUICK LOAD® SUTURE 2-0 PDO, Monofilament, Absorbable, 53" | Box of 12 Sutures (1 Suture per Pouch) | |

KNOTTING

| TK® PRODUCT ORDERING | | SUPPLIED: STERILE | | INDIVIDUALLY PACKAGED |
|--|------------|-----------------------------|-----------------------------|-----------------------|
| | REORDER | PRODUCT | DESCRIPTION | |
|  x 6 | REF 030404 | TK® Ti-KNOT® DEVICE | Box of 6 DEVICES | |
|  x 12 | REF 030510 | TK® QUICK LOAD® UNIT | Box of 12 QUICK LOAD® UNITS | |

5:00pm PN08XXXXA DRAFT Confidential ARR 05.23.23



[Scan to Learn More](#)

LSI SOLUTIONS®

7796 VICTOR-MENDON RD | VICTOR, NY 14564 USA | CUSTOMER SERVICE: 866.575.3493 | TECH SUPPORT: 866.428.9092 | FAX: 585.742.8086 | WWW.LSISOLUTIONS.COM

The LSI logo, LSI SOLUTIONS, JNW UrTrac, Monoglide, RD, RD180, Strongsorb, TK, Ti-Knot, and Quick Load are trademarks and registered trademarks of LSI SOLUTIONS, Inc. Copyright © 2023, LSI SOLUTIONS®. All Rights Reserved.